

Carepeutic®

Shiatsu Foot Massager

Instruction Manual
KH385L-B




OPERATION INSTRUCTIONS

TO START

- Connect the massager to power socket. When you plug in power cord, the On/Off LED light is on and the kneading board moves for about 5 seconds, this is to initiate the control system. After 5 seconds, the On/Off LED light is off.
- Press the On/Off button to turn on the massager. The massager is now ready to use.
- You can now choose auto program or customized program for massaging. Its default message program is P1.
- NOTE: For the first time use, you may hear gear sound from the message board during massaging. But the sound will disappear after a while.

AUTO PROGRAM

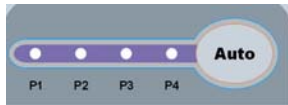

- Simply press the Auto button to choose programs, from P1 to P4. Each program lasts for 15 minutes. You'll hear an alert sound when you choose different programs (from P1 to P4).
- Its preset message time is 15 minutes for each program you choose (P1 to P4).
- If the massager is idol for 5 minutes with power on, it'll shut down automatically. Idol means that you turn on the power by pressing the On/Off button, but you do not choose any program from control panel.

General Relaxation

The following three programs concentrate on the areas of the foot that correspond to the following reflex point in the body:

- P1: Eyes, ears, neck, shoulders, spine.
- P2: Stomach, intestines, colon.
- P3: Legs, lower abdomen.
- Note: you can rarely feel the difference between P1 to P4 because it simply adjusts its strength on its reflex points only.

CUSTOMIZED PROGRAM



- Custom function is for you to choose the specific area that you want to massage on your feet (toe, arch, sole).
- Its preset message time is 15 minutes for each area you choose.
- You'll hear an alert sound when you choose different areas (toe, arch, sole).
- Note: you can rarely feel the difference between these 3 programs because it simply adjusts its kneading/rolling pattern on its reflex points only.

MANU PROGRAM



- The Manu button is just like a switch for you to turn on the Speed and Direction function.
- The "Speed" and "Direction" function can be used only after you choose "Manu" button.
- After you press the Manu button, you can now press Speed or Direction to choose massage speed intensity and or directions.
- When you press the Speed button, you can hear the "bee" sound, the Speed is the highest or lowest when you cannot hear the "bee" sound.
- Changing the direction is only for its moving kneading board pattern. It is not for clockwise or anti-clockwise direction.

TURNING OFF THE MASSAGER

- When the massager is turned on with Auto Program, it will turn off automatically after 5 minutes if it is not in use (if you do not press other buttons after finishing one program).
- If you press the On/Off button during massaging, the massager will be in pause mode. Simply press the On/Off button again and the massager will continue to do massaging.
- If you do not press any other buttons during pause, the massager will shut off automatically in 10 minutes.
- When you want to clean the fabric cover or the massager, you must unplug the power cord from power socket first.
- If you want to turn off the massager any time, press and hold the On/Off button for 5 seconds, and the LED light and power will be off completely. (After the LED light is off, the kneading board may continue to move for 2 to 3 seconds so that the kneading board can move back to its default position).

CAUTION (care and maintenance)

- When using electrical products, specifically when children are present, basic safety precaution should always be followed, read the instruction manual thoroughly before using this product.
- The appliance should never be left unattended when it is connected to power source. Unplug the power cord from power outlet when it is not in use.
- Always unplug the appliance's power cord from the electrical outlet immediately when it is not in use and or before cleaning it.
- Close supervision is necessary when this appliance is used by, close or next to children, invalids, or disabled person.
- Do not operate this appliance if it has a damaged power cord or electric plug; if it is not working properly; and if it has been dropped or damaged or immersed into

-
- water. Return the product to service center for examination and repair. Never try to repair this device by yourselves.
- Keep power cord away from heated or wet surfaces.
 - Do not operate this device under blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to person.
 - Do not carry this appliance with the attached power cord or use the power cord as a carrying handle.
 - This device is NOT for outdoor use.
 - Do not spray water into the massager or the fabric covers. Please remove the fabric cover for cleaning.
 - To disconnect the power source, turn off power from On/Off button and then remove the plug from power outlet.
 - Do not attempt to stand up on the massage board. Use it only when you are being seated on a chair or sofa.
 - Connect this appliance to proper electric outlet only (110/120V).
 - Do not use this massager with wet feet.
 - Do not use abrasive cleaning agents such as benzene, thinner or other solvents to clean the massager. Use only dry cotton cloth to clean the surface of massager.
 - Do not put this massager next to heater or stove or direct sunlight.

Special Warning

- Do not use the device in the bathing room or where you take show. Use it only in a dry area. Keep it away from humid area, wet floor. Do not spray water onto the control panel or the outside case areas around the appliance.
- This appliance is not intended for use by person (including children) with reduces physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction in regards to the use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Do not use this device continuously for more than 45 minutes at a time. Extensive use could cause overheated on the device. Should this occur, turn off the power completely and allow the unit to cool down before using it again.
- Never use the appliance while sleeping or drowsy. Avoid smoking or drinking during massaging.
- This product should NEVER be used by any individual suffering from any ailments that would limit to altogether inhibit the user's capacity to feel or have any sensation in the body part that the product is being used for.
- Consult your doctor before using this unit if you have diabetes, as diabetes can contribute to foot problems including decreased feeling or circulation to the feet. Given the existing risk of foot problems for diabetics, this unit must not be used prior to consulting your doctor.
- This appliance is not intended for use by person (including children) with reduces physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction in regards to the use of the appliance by a person responsible for their safety. Children should be supervised

to ensure that they do not play with the appliance.

Warning: Any individual who is pregnant, has a pacemaker, suffers from diabetes, or suffers from malignant tumors, excessive bleeding or tuberculosis should consult with a physician before using a massaging device designed for home use.

SPECIFICATIONS

Power voltage: 110/120V

Power consumption: 40W

Auto Massage timer: 15 minutes

Speed: High, Medium, Low

Dimension: 65cmx39cmx31cm

Weight: 9.5KGs

Storing temperature: 15C to 85C

Normal operating Temperature Range: +10c to +40c



Note: This product is not intended to diagnose, treat, cure, or prevent any disease.

Limited Warranty:

Please contact seller for warranty. Limited warranty from manufacturer only for original purchaser against manufacturing or workmanship defect under normal use. Due to the hygiene or rehabilitation nature of this product, it is non-returnable.

Carepeutic® is a registered trademark, all rights reserved, © Q.L.M.

www.ideatomorrow.com

Powered by Carepeutic® Global Lab