

SHARPER IMAGE®

LED HOME ELECTROLYSIS KIT

Item No. 201622

Quick Start Guide




Thank you for choosing the LED Home Electrolysis Kit from Sharper Image. Please read this Quick Start Guide to get a general idea of how this product works. For additional instructions and important safety information, be sure to read the entire user's manual in the package.

FEATURES

- Painless hair removal at home
- Safe for all skin tones and all hair types
- Cordless and portable
- Includes small roller tip for facial hair and eyebrows
- Includes large roller tip for upper lip, bikini line and other body hair
- Includes 3 self-adhesive gel pads (small, large and upper lip) for treating larger areas
- Includes treatment tweezers for removing single hairs
- Includes conductive gel, cleanser and carry case
- Operates on 2 included AAA batteries

IDENTIFICATION OF PARTS

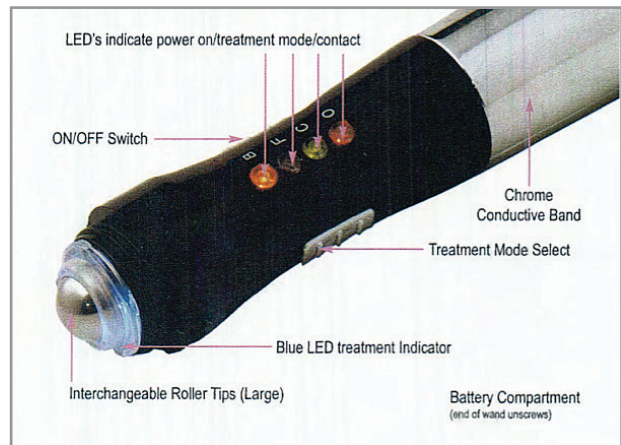

BATTERY INSTALLATION

- Install the 2 AAA batteries (included) by twisting the battery cap counter-clockwise.
- Insert the batteries according to the correct polarity ["+" or "-"].
- Replace the cap.

REMOVING FACIAL HAIR WITH INTERCHANGEABLE ROLLER TIPS

Select one of the two included roller tips. The SMALL roller tip is for facial hair and unwanted eyebrow hair. The LARGE roller tip is for the upper lip. To change tips, hold the clear plastic ring and gently pull toward you. Then insert the other tip into the socket and press gently into place. [The device comes with the large roller tip already installed.]

1. Remove makeup and clean the area thoroughly with the included pre-epilation cleanser. Treatment is most effective after a bath or shower.

REMOVING FACIAL HAIR WITH INTERCHANGEABLE ROLLER TIPS [CONT.]

2. Apply a small amount of the included conductive gel to area you will treat. This area should be no more than 0.3-0.5 sq. in. [2-3 sq. cm].
3. Slide the Power switch to ON.
4. Slide the Treatment Mode switch until "F" [for Face Mode] is illuminated with a YELLOW LED.



5. Holding the device in your hand like a pen, place the roller tip on the area to be treated. Make sure your finger and thumb or palm are touching the

chrome conductive band on the device. Gently roll the roller tip back and forth over the area, making sure the roller tip does not lose contact with the skin. One beep will sound, indicating that proper contact has been made. The GREEN contact LEDs will light up, and the BLUE LEDs around the roller tip will light up during treatment. [The GREEN and BLUE LEDs will flash while the device is in rest mode.]

The treatment cycle is as follows:
20 seconds on >> 10 seconds rest >>
20 seconds on >> 10 seconds rest >>
20 seconds on >> 10 seconds rest >>
3 beeps >> TURNS OFF AUTOMATICALLY

NOTE: If you feel a slight tingling sensation, or the GREEN LED does not light up, cleanse the treatment area or add more conductive gel.

6. Gently tweeze hairs in the direction of growth.
7. After treating one area, move to another area, following the steps above.

8. Apply your normal moisturizer after treatment, as electrolysis can dry your skin.
9. Remove the batteries and clean the roller ball with a tissue or cotton ball. Wipe the wand with a tissue or soft cloth. **DO NOT IMMERS THE DEVICE OR THE ROLLER TIPS IN WATER.**

NOTE: For FACE MODE, the recommended DAILY treatment is two 90-second cycles per area, per day. Do not exceed this daily amount. The recommended WEEKLY treatment is two 90-second cycles per area, per day, for up to 5 DAYS in a given week. Do not exceed the weekly maximum of 10 cycles per week.

This versatile device also removes body hair with the self-adhesive gel pads and single hairs with the included treatment tweezers. Consult your manual for additional instructions and important safety warnings.

Remember, hair removal is a cumulative process. The longer you use this product, the more effective it will be. Hair grows in a cycle of 8-12 weeks. With regular use, you should see results in about 60 days.

CUSTOMER SERVICE

If you require any assistance with your LED Home Electrolysis Kit, please call 1 (877) 210-3449.

SHARPER IMAGE®