



Personal Massager

PL021



Operating Manual

Introduction

Prospera's PL021 Personal Massager, a user friendly product, can relax muscles and relieve pain on your waist, shoulders, necks, arms, legs and abdomen with its powerful and high-speed vibration massage. You can adjust the massage intensity of PL021, which includes four operating modes and a soft & waterproof head. It is a perfect personal massager for use at home and in office.

Features of PL021 Personal Massager:

1. Four speeds with different intensities
2. Soft and waterproof head for gentle massage
3. Muscle relaxation and pain relief on various body areas

Content

CONTENTS

Safety Warnings.....	3~4
Set up.....	5
Operating the product.....	6
Trouble shooting.....	7
Consumer warranty.....	8
Positions.....	9
Warranty card.....	10

Safety Warnings

When the Personal Massager is not in use, always unplug it from the socket. This action prevents electric shock and damage to the product.

Warnings:

Consult your doctor prior to using the product if you have heart, skin and/or other health-related issues.

Please take the following precautions to prevent electric shock, burns, fires and other injuries.


1. Do not leave this product unattended, and unplug from the socket when it is not in use.
2. Do not operate this product in an excessively humid environment, and avoid direct sunshine.
3. This product is intended for adult only, and do not use the product outdoor.
4. Keep the power cord far away from heat sources, and do not disassemble the product unless authorized.

Attention:

1. Do not let babies and children play with the product.
2. Start massaging from the lowest intensity of product. If you feel uncomfortable, reduce the intensity or stop using the product immediately.
3. Recommended daily usage is about 20 minutes.
4. Do not use it while you are asleep.

5. Do not use the product on eyes and/or any sensitive areas.
6. Do not use the product under the blankets and the pillow to avoid overheat.
7. Do not use the product when the skin is swelling, inflamed, and/or abnormal.

Set up

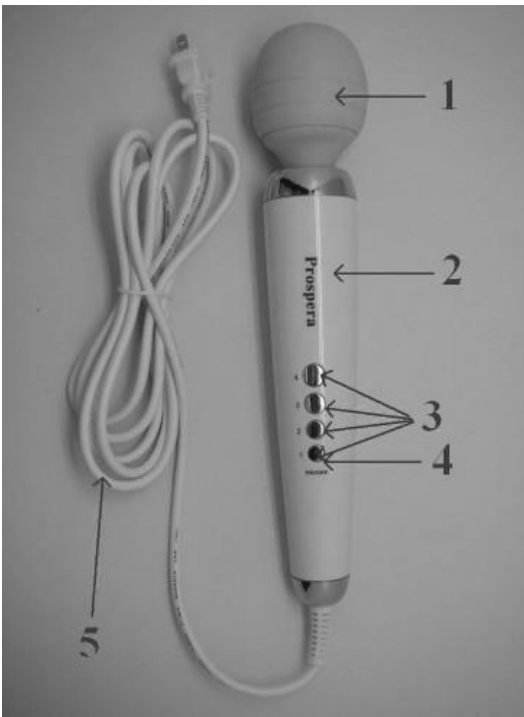

1. vibration head
- 2 .main body
3. speed buttons
4. speed indicator
5. Power cord

Product Specifications:

Unit size	12.2"×2.4"×2.4"
Unit Weight	1.44lbs
Length of the cable	98.5"
Alternating Voltage	110V
Frequency	60HZ
Vibration level and speed	Level 1: 3500 rpm Level 2: 4200 rpm Level 3: 5000 rpm Level 4: 6800 rpm

Operating the product

1. Press the ON/OFF button to turn on the massager, which is at the lowest speed (Level 1)
2. Press Buttons 2-4 to adjust the vibration speed of the massager (Level 4 is the strongest level)
 - Level 1 and switch on/off
 - Level 2
 - Level 3
 - Level 4
3. After the level of intensity and speech is selected, put the massager on the body areas, where massage is needed.
4. The massager will turn off automatically after 20 minutes of consecutive usage
5. Press the ON/OFF button to turn off the massage
6. Recommended time of usage: 20 minutes

Troubleshooting

Prospera's PL021 Personal Massager is a user friendly product. If you follow the operating manual, this product will last for a long time.

However, if the following situations happen, please take the recommended troubleshooting procedures:

1. If the soft head of the product becomes loose, stop using it immediately, and report to the dealer or Prospera.
2. If the powder cord is worn off or the product does not work properly, unplug the product from the socket. Contact Prospera if it is under warranty.

Consumer Warranty

Prospera offers a 90-day standard warranty to the massage products. Labor and parts are included during warranty period. Consumers are required to provide model number and dated proof of purchase (sales receipt) when contacting Prospera for the warranty service.

Non-warranty service is provided on a “per incident” basis. The consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and “per incident” fees will be charged for non-warranty repairs and support calls.

Positions



Warranty Card



Warranty Registration Form

Product Name: _____

Product Model: _____ Purchase Price: _____

Date of Purchase: _____ Retailer: _____

First Name: _____ Last Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Phone: _____ Email: _____

Signature and date: _____

Please mail to:

Prospera Corporation
405 Boulder Ct, Ste 500
Pleasanton, CA 94566

Telephone: 925-225-0888
Fax: 925-225-0660
Email: sales@prosperacorp.com

