

# SHARPER IMAGE®

## 3-IN-1 NECK THERAPY WITH REMOTE

Item No. 205388

Owner's Guide



Thank you for purchasing the Sharper Image 3-in-1 Neck Therapy with Remote. Please take a moment to read this Owner's Guide and store it for future reference.

### FEATURES

- Soothes and relaxes the neck
- Combines low-frequency electrical impulses, Far Infrared Ray (FIR) heat, vibration massage and embedded magnet therapy
- Far Infrared Ray (FIR) heat penetrates muscle and bone tissue without getting hot to the touch
- Manual and automatic massage modes (6 preset modes)
- 16 electrical intensity levels
- 3 vibration intensity levels
- 3 heat intensity levels (High/Low/Off)
- 6 electrical impulse modes
- 6 vibration modes
- Automatic shutoff timer (5-30 minutes)
- Comfortable, ergonomic design
- Cordless design for portability
- Wireless remote control with backlit LCD screen



### CHARGING THE 3-IN-1 NECK THERAPY WITH REMOTE

*IMPORTANT: Please charge this device before using it for the first time.*

- Plug the adapter into the DC power supply jack on the device.
- Plug the other end of the adapter into an AC outlet.
- The LED light will flash quickly while it is charging.
- The LED light will flash slowly once it is fully charged.





*NOTE: A full charge takes approximately 2 hours.*

## USING THE 3-IN-1 NECK THERAPY WITH REMOTE

- Wipe away any perspiration from your neck.
- Moisten the neck ring with a slightly damp towel.
- Place the device around your neck [with the opening in front].
- Set the function switch to the ON position.
- Press the Power button  to begin operation. *NOTE: The LCD screen on the remote displays the status as it enters the default mode [automatic mode]. See description of modes below.*
- Press the Power button  once more to cease operation.

## USING THE WIRELESS REMOTE



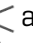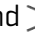


The wireless remote control uses infrared light emitted from the bottom of the controller. The infrared light must make contact with the infrared receiving window on the front of the device. See illustration below.

- Press the  and  buttons to scroll up and down.
- Press the  and  buttons to decrease or increase values.

*NOTE: If the wireless remote stops working, replace the battery. Use a fresh CR2032 button cell battery.*



## MESSAGE MODES

- Press the  or  buttons to enter the mode menu . Then, press  or  to select one of 6 message modes:

MODE 1 [Automatic mode]: Vibration + Heat / Heat + Electrical Impulse / Heat only / Vibration only / Electric Impulse only

MODE 2: Vibration + Heat

MODE 3: Heat + Electrical Impulse




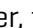

MODE 4: Heat only

MODE 5: Vibration only






MODE 6: Electrical Impulse only

*NOTE: Do not change modes while this device is in use. Turn it off first and then restart it.*


## TIMER SETTINGS

- Press the  or  buttons to enter the timer settings . Then, press  or  to set the automatic shutoff timer, from 5 to 30 minutes.


## ELECTRICAL IMPULSE SETTINGS

- Press the  or  buttons to enter the electrical impulse settings . Then, press  or  to select one of 6 different impulse modes.

## SELECTING INTENSITY LEVEL

- Press the  $\wedge$  or  $\vee$  buttons to enter the intensity settings . Then, press  $<$  or  $>$  to select the intensity of the electrical impulses [16 levels] or vibrations [3 levels].



## VIBRATION MODES

- Press the  $\wedge$  or  $\vee$  buttons to enter the vibration mode settings . Then, press  $<$  or  $>$  to select one of 6 different vibration modes: Auto / Vibration / Kneading / Beating / Press / Rolling.

## SELECTING HEAT LEVEL


- Press the  $\wedge$  or  $\vee$  buttons to the heat settings . Then, press  $<$  or  $>$  adjust the heat to High [H], Low [L] or Off [0].

## DESCRIPTION OF PARTS

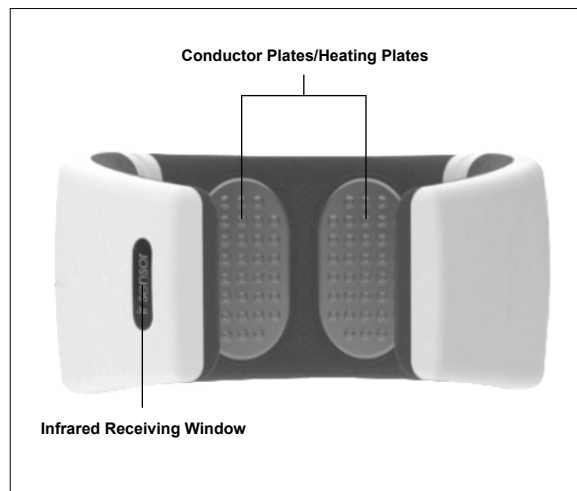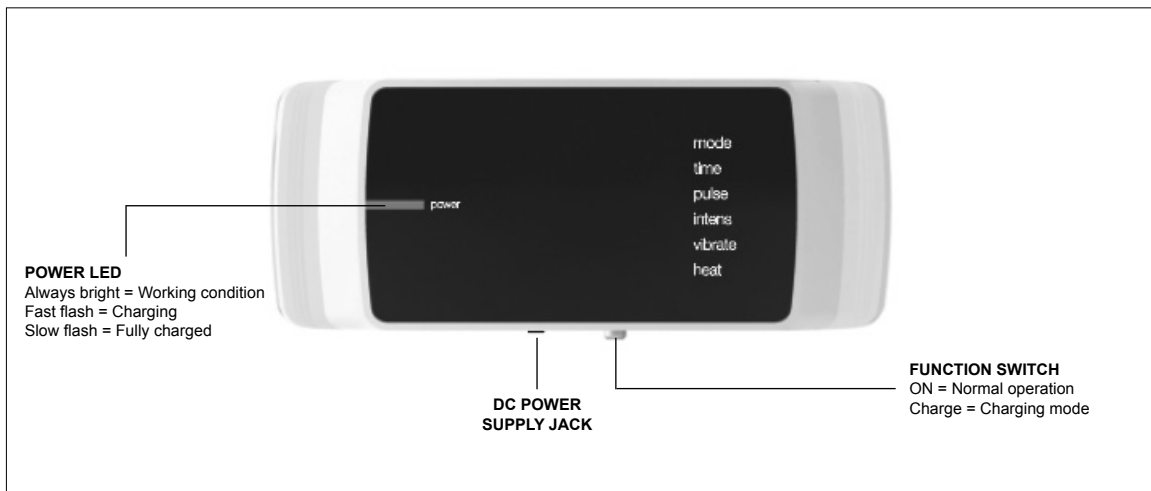

## DESCRIPTION OF REMOTE

- A. LCD Screen
- B. Mode
- C. Timer
- D. Electrical Impulse
- E. Intensity
- F. Vibration
- G. Heat
- H. Scroll Up
- I. Scroll Down
- J. Decrease Value
- K. Increase Value
- L. Power Button
- M. Infrared Signal Emitter
- N. Battery Compartment
- O. AC Adapter

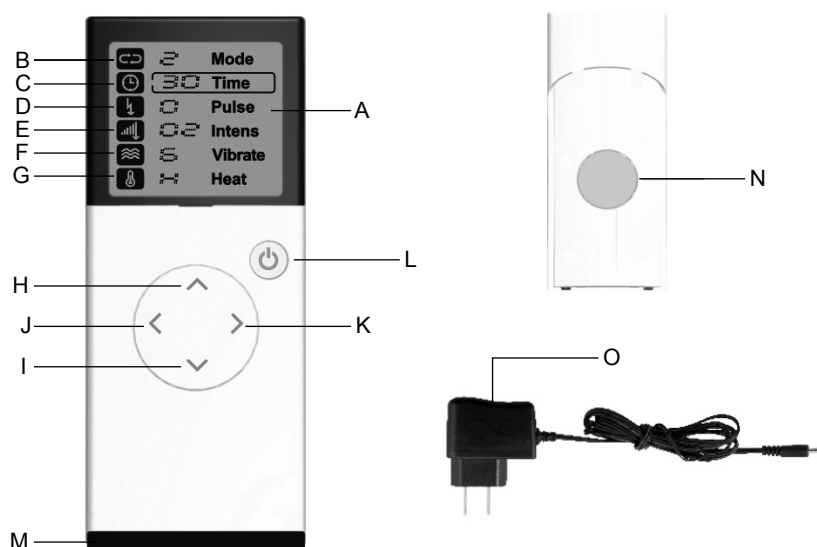

## TECHNICAL SPECIFICATIONS

Model:	URG-13178
Weight:	630 g
Power [device]:	6V DC ≤ 6W
Power [remote control]:	CR2032 button battery 3V DC
Electric impulse frequency:	1-1000Hz
Width of electric impulse:	20-400us
Wave of electric impulse:	square wave

## WARNINGS

- Do not use this device if you have a pacemaker.
- Do not use this device if you have sensitive skin or a heat sensitive condition.
- Do not change modes while this device is in use. Turn it off first and then restart it.
- Pregnant women and cardiovascular patients should consult their doctor before using this device or any at-home massage device.
- Do not use this device while sleeping.
- Do not use this device in the bath or shower.
- This device is not a toy. Keep out of reach of children at all times.
- An increased intensity of pressure or power does not result in a more beneficial therapeutic effect.
- Do not open or dismantle this device. There are no user-serviceable parts.
- Do not immerse this device in water or any other liquid.

## WARRANTY / CUSTOMER SERVICE

Sharper Image branded items purchased from SharperImage.com include a 1-year limited replacement warranty. For Customer Service, please call 1 (877) 210-3449.

## SHARPER IMAGE®

The SHARPER IMAGE® name and logo are registered trademarks.  
©Sharper Image. All Rights Reserved.