

Operation Tips:

- Noise after pressing "Pressure" button: after pressing the "Pressure" button, you will hear the gear sound--that is the nature of modes switching, then press the "Inflatable" button to adjust the pressure until it reaches the pressure level you like, and the gear sound will be gone.
- Noise after choosing the Slow Speed Mode: when you choose the Slow Speed mode, it is the nature to hear the sound from the 192 acupuncture nodes and that sound will disappear after switching it to Fast Speed mode.
- Please do NOT stand up on the 192 nodes massage board. It'll cause damage to the panel inside.
- No power: Please turn on the Power Switch Button located on the Bottom of the Massager before pressing the On/Off Button from the touch panel of the massager.



Acu-Point Therapy Foot Massager

The First Foot Massager Combined with Chinese traditional Acupuncture Therapy and Heated Air Massaging to Improve Blood Circulation, relief painful, sore, tired feet.



EMK-730C

INTRODUCTION

This foot massager applies acupuncture technology with acupuncture nodes combined with air pressure massaging and infrared heated air to improve blood circulation.

FEATURE

- Adjustable 4 Unique airbags strengthen finger pressing massage.
- Preset airbags pressure in 3 levels of intensity.
- Manual inflatable and venting adjust air pressing.
- 192 acupuncture nodes create reflexology.
- Adjustable finger pressing speeds in 3 levels.
- 24 infrared lights enhances blood circulation
- Time settings are 10, 15 (auto) and 20 minutes




192 Acupuncture Nodes



4 Airbags



24 Infrared Lights

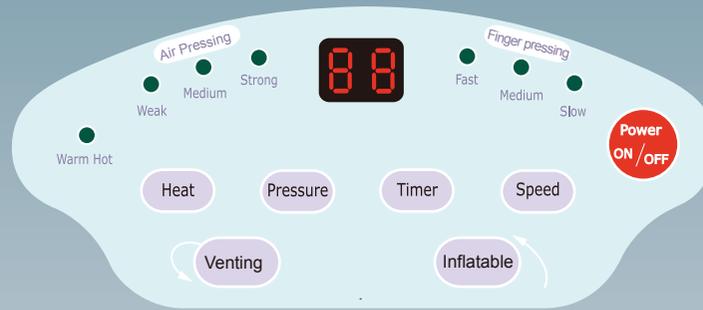

NOTICE

Please turn on the Power Switch located on the bottom of the foot massager after plugging in the power cord.

Please turn off the Power Switch located on the bottom of the foot massager after plugging in the power cord or when it is not in use for long time.

OPERATION

1. After you turn on the Power Switch located on the bottom of the massager, you'll hear the "beep" sound once.
2. Power On/Off
Press the "On/Off" button once to initiate the device and it will get into standby mode with automatic preset timer of 15 minutes showed on the LED panel.
3. Speed
Press "Speed" button to choose 3 different speeds (Slow/Medium/Fast) showed on LED panel.

4. Timer
Press "Timer" button to choose 3 presets operation timers (10/15/20) minutes showed on LED panel.
5. Pressure
Press "Pressure" button to choose 3 preset levels of air pressing (Weak/Medium/Strong) showed on LED panel.

Technical Tips: after press the "Pressure" button, you will hear the gear sound, then press the "Inflatable" button to adjust the pressure, and the gear sound will be gone.

6. Inflatable
Press "Inflatable" button repeatedly until reaching the pressure level you want.
7. Venting
Press "Venting" button repeatedly to release air pressure until reaching your desired level.
8. Heat
Press "Heat" button to turn on or turn off the heat function showed on LED panel.
9. Press the "On/Off" button to turn on or turn off the device.
Turn off the Power Switch located on the bottom of foot massager to stop operation.

Specifications Voltage:
AC100-240V/50-60Hz
DC24V 2 A
Power: 48W
Weight: 6Kgs