


operation instructions

Five Elements Meridian Brush

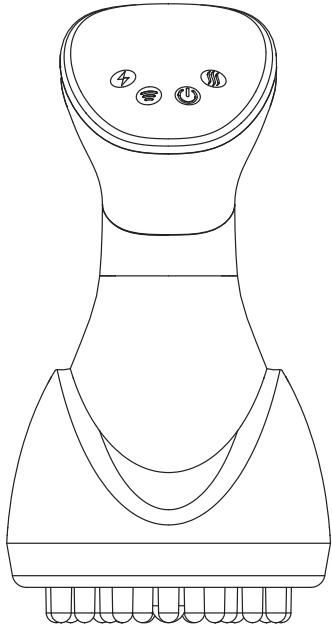

Note: This product cannot be used as a substitute for drugs or medical devices

The pictures are for reference only, and the product is subject to the actual product!
Please read the instruction manual carefully before use


preface

To Dear Customer:

Thank you very much for purchasing our company's Five Elements Meridian Brush. To ensure your correct and safe use, Please read the user manual carefully when using this product. Lend or transfer to others for use, be sure to keep this. Please forward the instruction manual together.

This product is not a medical device and is not recommended for medical use. The product uses electromagnetic pulse (TENS) technology to simulate massage techniques to relax muscles, and any electric flu during the process is normal

Introduction to Operation Panel

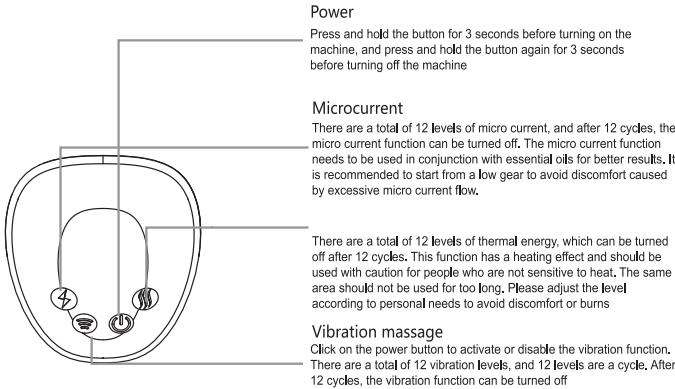

When charging, the indicator light is red and the fully charged indicator light is green. When charging, please use 5V/1A appliances and keep away from flammable and explosive materials such as wood, cotton, paper, etc. Incorrect use of chargers may lead to product malfunctions.

Reminder: The instrument has been set to automatically shut down on a scheduled basis for safety reasons, in order to prevent potential hazards of prolonged heating caused by transportation or accidental operation.

Attention: The product cannot be turned on and used during charging. Please unplug the power to turn it on. When not in use for a long time, please store it fully charged and charge it once a month, otherwise it will damage the battery and make it unusable or shorten its service life

note

- Do not use voltage other than DC5V for charging, and do not use it during charging.
- Do not use in hot places such as bathrooms, do not immerse yourself in water or get close to water sources to avoid short circuits or electric shock.
- When cleaning, in case of malfunctions during use, or during use, the machine must be turned off and the power cord must be unplugged. 4. In thunderstorm weather, please do not connect the power to charge to avoid excessive current, damage to the massager, and even personal injury.
- Do not charge for a long time. After fully charging, unplug the power in a timely manner to avoid damage or safety issues caused by overheating of the battery.
- Please use neutral cleaning to wipe the dial. Do not use automotive repair, oil refining, or highly acidic or alkaline corrosive cleaning agents for cleaning.
- After use, place the massager in a safe, dry, and cool place.
- Please avoid collision or heavy pressure, and avoid hard and sharp objects scratching the surface of the machine to avoid unnecessary malfunctions.
- When not in use for a long time, please store it fully charged and charge it once a month, otherwise it will damage the battery and make it unusable or shorten its service life.

Environmental information

Environmental protection: The built-in rechargeable battery contains substances that may pollute the environment. Before discarding the meridian drama, be sure to remove the battery or send the rechargeable meridian brush to the official designated recycling station for disposal. You can also send the discarded charging meridian brush to the company's designated maintenance department, which will decompose and recycle the battery.

Disassembling the battery: Cut off the power and let the charging meridian brush run until the battery runs out. Pry off the black panel from the button position, remove the screws, and then remove the battery.

 **Only remove the battery when discarding the charging meridian brush.**

Kind reminder:

This product should not be disposed of together with household waste. It should be concentrated in appropriate institutions for recycling and recycling. The lithium-ion batteries used in this product are valuable resources that can be recycled and reused. After use, please hand over this product to professional personnel to remove the lithium-ion batteries and hand them over to professional institutions that can recycle and regenerate rechargeable batteries for disposal.

Product parameters

Product Name: Five Elements Meridian Brush

Rated voltage: 5V ---

Rated power: 12W

Execution standards: GB4706.1-2005 GB4706.10-2008

Sole: Massage from heel to toe

Lower limbs: undocking from bottom to top (i.e. massaging towards the internal organs);

Upper limbs: unblock from top to bottom;

be used:

Abdomen: The abdomen can be cleared from top to bottom (for those with visceral prolapse, it should be cleared from bottom to top), or spiral therapy can

the body.

Chest: Use a heavy brush to clear the midline of the chest from top to bottom. The chest is bounded by the Ren meridian on both sides of the anterior midline of

Back: Smooth from top to bottom

Neck: Brush from top to bottom, you can often scrape and wipe the neck;

Head: Use a meridian brush from the front of the forehead to the back of the

the limbs:
times, and scrape until the head feels hot;
brain, without the need to apply lubricant randomly. Clear each area about 30

Usage: From top to bottom, massage the left side of the body first, then the right side of the body, starting from the head, neck, back, waist or neck, and back to

Massage methods

Tips

- When cleaning the instrument after use, it should not be washed with water. You can use a pure cotton cloth with clean water and wipe it with household alcohol for future use.
- It is recommended to use each part for about five minutes each time, and a single use time of about ten minutes is recommended.
- There are 12 adjustable gears for the left micro electric and right thermal energy, and it is generally advisable to turn them to the left micro electric and right thermal energy, and use time of about ten minutes is recommended.

warn

Disalbed in the following situations:
1. Patients with skin diseases and skin allergies
2. Pregnant women and individuals in their menstrual period
3. Heart disease patients
4. Patients with abnormal blood pressure
5. Osteoporosis patients
6. Malignant tumor patients
7. Cardiovascular patients
8. Patients with spinal deformities or abnormalities
9. Individuals with foreign body implants in their bodies
10. Patients with neurological abnormalities in the brain
11. Patients with acute suppurative inflammation
12. Patients with bleeding tendencies

Please consult a doctor before use in the following situations, otherwise accidents or physical discomfort may occur :
1. Patients with diabetes
2. Patients who cannot express their wishes
3. Patients who are not sensitive to heat
4. Patients receiving medical treatment
5. People who feel unwell
6. Patients with a body temperature exceeding 38 degrees Celsius (fever period)
7. Patients with infectious diseases
8. Patients taking medication
9. Beneficiaries with metal allergies
10. Patients with skin perception disorders
11. Patients with skin abnormalities

Safety precautions

4. Before using the instrument, apply a hot towel to the area that needs to be massaged, which can moisturize the skin. Because dry and irritable skin can cause uneven penetration, resulting in a stinging sensation.

1. If there is varicose veins in the body, it should not be used.

2. During pregnancy or menstruation, the use should be temporarily suspended.

3. Do not use during charging and keep away from flammable materials.

4. Do not use for a long time, the best massage time is 15 minutes.

5. Children are strictly prohibited from using it, and special users must use it under adult supervision.

6. It is strictly prohibited to dismantle the machine without authorization, otherwise any malfunction or danger will occur, which is not covered by the warranty. Our company will not be responsible for any adverse consequences caused by this.

7. When using, it should be brushed in one direction, not one up and one down, to heat must pay attention when using it.

8. The surface of the appliance generates heat, and people who are not sensitive to heat must pay attention when using it.

9. When massaging, the micro electric brush head should not massage near the heart, otherwise it will increase the risk of heart fibrillation.

10. The product is recommended for one customer to use. When using it for cross users, please clean the brush head with disinfectant before use.

11. This product contains batteries. Please do not disassemble, modify, flood the product, or throw it into fire, as this may cause damage or explosion to the product.

12. If you experience any discomfort during or after use, please stop using immediately and consult a doctor.

13. During massage, the micro electric brush head should avoid contact with wounds and scars.

quality assurance

We hereby guarantee that the product has excellent quality, reliable performance, and is safe and practical. To address customer concerns, our company implements quality responsibility and after-sales service, which can be processed with this card.

1. If the main components of the product malfunction during normal use within one year, free repair and replacement of the components will be provided.

2. When repairing, payment vouchers and guarantee cards should be brought to determine the warranty period.

3. Product warranty period of one year

take care

1. Abnormal use and human damage are not covered by the warranty.

2. If the customer disassembles the machine on their own, the warranty qualification will be cancelled.

Purchase date: _____

Purchase location: _____

certificate qualified

Inspection results: _____
examination clerk: _____
Inspection date: _____

Special declaration

1. If there are any technical changes to this product in the future, no further notice will be given

2. If there is any change in the appearance or color of the product, the actual product shall prevail