
# BCE1说明书 105g双铜,风琴折,90\*150mm,8P,单黑色印刷,过光油

# LUNA'MOUR

# CelluFade BodyBliss

EMS + Vacuum Sculpting + Red Light Therapy Body & Face Device

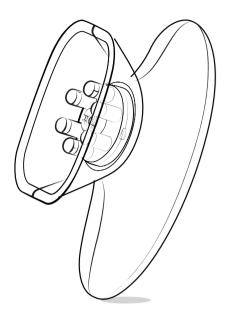

Instruction Manual Model: LL-CELLU


Please read this manual carefully before usage.

# MEET THE CELLUFADE

Meet the CelluFade BodyBliss device! Using Vacuum Sculpting Therapy, Electrical Muscle Stimulation, and red light therapy, the CelluFade is designed to target lumps and bumps by helping break down fat cells, stimulate collagen production, and reduce the appearance of cellulite.

With 3 technology modes and 3 intensity levels, the CelluFade BodyBliss offers a customized treatment session that can be adjusted to your comfort level. This 3-in-1 device can be used on stomach, legs, buttocks, arms, and face.

The CelluFade is easy to use and will become a staple in your skin care routine.

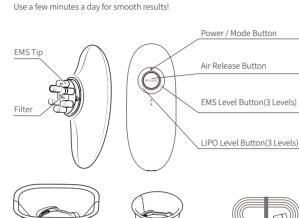

Small Suction Cup Type-C Charging Cable \*EMS: Electrical pulse current that stimulates muscle contractions

\*LIPO: Stimulating muscle movement for a lifting and shaping effect.

## **PRECAUTIONS**

Please read the instructions carefully before using the device.

1 Keen the device away from water

2. Do not immerse or drop into the water or any other liquid.

3. If the device falls into water, unplug it immediately. Do not place your hand in

4. Use the device only for its intended purpose, as indicated in this manual. Do not use any accessory which is not recommended by the manufacturer 5. Do not drop the device or insert any object in the openings.

6. Do not overuse, maximum 5 minutes on treatment area. Do not stay on one area for more than 5-8 seconds.

7 For people with loose muscles, or skin this device may leave a mark after. usage. In case bruising occurs, stop usage immediately. Bruising should disappear gradually.

8. If you have a medical condition, or are unsure if this device is right for you, please consult your doctor before using the device.

9. Do not open casing of the device.

10. Never leave the unit unattended, especially when children are present 11. The device should only be used with the power supply included with the

## **HOW TO USE**

### Charging

1. When battery is running low to last one cycle (10 minutes usage), all LED lights will flash 5 times as a warning, then, enter into regular working mode for

2. Plug the Type-C end of the cable into the charging port on the device.

3. Plug the USB end of the cable into a standard USB charging plug. 4. During charging, the LEDs on main unit will illuminate in sequence, and form

a full circle when fully charged.


5. A full charge requires up to 2  $\sim$  4 hours and will provide up to 60 minutes of usage. Once fully charged, unplug the adaptor from the socket.

### Instructions for usage

1. Choose the suction cup according to the size of the area to be treated.We recommend using the LARGE CUP ideal for ABDOMEN, WAIST, BACK, HIPS, LEGS, BUTTOCKS. Use the SMALL CUP ideal for ARMS and CALVES.

2. Gently place the suction cup on the main unit at a 30 degree angle. Remember, the "pointing mark" on cup is to the lower right corne

3. Twist the cup clockwise, until you hear a "click" sound, "Pointing mark" is in the middle of the device.

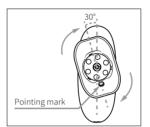



## **HOW TO USE**

4. Press Power button to start the device. Press once again to switch mode between Body Mode and Face mode.

5. Under Body Mode, press LIPO button to adjust suction level, lights will turn on accordingly. The suction level will change in the following sequence by further pressing the LIPO button. Initially, start at level 1 and gradually increase suction intensity according to your comfort.


Level 1  $\rightarrow$  Level 2  $\rightarrow$  Level 3  $\rightarrow$  Suction OFF

6. Press EMS button to adjust EMS level with lights indication illuminated. The EMS level will change in the following sequence by further pressing the EMS button. Initially, start at level 1 and gradually increase level as the body becomes accustomed to EMS stimulation.

### Level 1 → Level 2 → Level 3 → Suction OFF

7. Under Face mode, only EMS is activated, the suction feature cannot be used Use the device without either cup

8. Once the settings are as desired, firmly place the cup on skin and massage the treatment area by keeping a slow and steady movement. Always massage in upwards direction, activating the air release button after each upward movement, before initiating another.

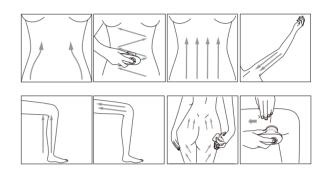



# **HOW TO USE**

9. Adjusting the suction speed level to your comfort is essential. When treating multiple body parts, start with the upper area and progress downwards. 10. Follow the recommended motions depicted in the accompanying images, and adhere to the suggested durations of 3-5 minutes per area, applying treatments every other day until desired results are achieved. Subsequently, reduce frequency to once every other week to maintain inch reduction RECOMMENDATION: It is ADVISED to apply a water-based gel a water-based gel such as aloe vera on the targeted treatment area to smoothly glide device over the skin and a stronger conductivity for EMS therapy.

11. CAUTION: Do NOT overuse, maximum 5 minutes in one area. Note: It is normal for skin to temporarily turn red after usage, which will disappear after 1 day of non-usage.

12. During usage, keep a slow and continuous movement. Simply hold the device firmly against your body to maintain the cupping suction. CAUTION: Do NOT stay in one position for more than 5-8 seconds. 13. When both LIPO and EMS are on, you should feel both suction of your skin and light electric pulses. If either becomes uncomfortable, press the RELEASE button for immediate stop. Adjust the level to a lower setting to re-start the

14. When both LIPO and EMS are on, you should feel both suction of your skin auto-shutoff after 10 minutes.

15. Twist cup anti-clockwise to remove from main unit.

# HOW TO USE

Cleaning

1. Make sure the unit is un-plugged and is powered off.

The filter in your CelluFade device should be replaced once a month. To replace filter, make sure device power is OFF. Remove old filter and insert new filter as shown in the image.



IMPORTANT: KEEP the unit away from water. Do not rinse or immerse the unit

with water.
IMPORTANT: Never use detergent to clean the main unit or suction, to avoid

After  $\overset{\smile}{\mathsf{usage}}$ , make sure to power off the unit and remove the adapter from the Store the unit away from areas with high humidity or high temperature.

This product should NOT be treated as household waste when you wish to dispose of it. Instead, it should be handed over to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring the device is disposed correctly, you will help to reduce negative consequences to the environment and human health, which may otherwise be caused by inappropriate disposal of the device. This recycling of materials will help to

Please include our standard disclaimer, warranty, and contact info as used in all

- 6 -

# Warnings:

# CONTRAINDICATIONS

No modification of this device is allowed

Ask your doctor prior to using if you:

Are taking prescription or over the counter drugs known to cause light sensitivity.

Are, or suspect you are, pregnant.

Have any suspicious or cancerous lesions.

Please read the instructions carefully before use and use this product in the

Please keep out of reach of children. This device is designed for use on adults over 18 years of age and not recommended for use on children under the age of 12.

Please do not use if you are experiencing the following: Heart disease, high fever, infectious diseases, malignant tumors, skin wounds, acute diseases, using medical appliances, skin diseases, pregnant women and/or breastfeeding, excessive sunburn, particularly sensitive skin, and those who are in

Please consult your doctor prior to use if you are, or suspect you are, pregnant. Please consult your doctor prior to use if you are, or suspect you are, pregnant. Please consult your doctor prior to use if you have, or suspect you have, cancer. Do not use device on open wounds. Do not use this device on skin that is painful to touch, raised (higher than the surrounding skin), crusting, inflamed, or bleeding. This may be a sign of a more serious condition and should be treated by your doctor.

- Do not place magnets, metal or conductive objects on the body of the device as it may cause malfunction.

- Do not rinse the instrument with water, use protective cloth provided to wipe clean after use.

- Do not wipe with alkaline detergent, alcohol, thinner, gasoline, nail polish

- If not in use for a long time, please make sure that the instrument is fully

- Please stop using it immediately when there is abnormal heat, so as not to cause burns and fires and other safety issues. - When the power cord is damaged, please do not use it to avoid short circuit and

- DO NOT submerge the device into water as it will malfunction and/or cause

- Do not place near flammable materials (spray, oil, natural gas, etc.)

- If product has abnormal noise, odor, or smoke, do not disassemble it at will, please contact the customer service email: hello@lunamour.co

-7-

# Important Safety Information:

- If the product is faulty, please contact the customer service center, please do not disassemble or modify the machine by yourself.

Do not operate device where it would come in contact with water. Do not use or charge device with power adapter that is damaged and/or faulty. Do not harshly bend or abuse the wire of USB charging cord. Do not attempt to take the device apart. Device is not serviceable by user. This will void any warranty.

> Stop using the device immediately if you think it is defective Please contact Client Care at: hello@lunamour.co

LUNA'MOUR is constantly improving its products with the latest technology. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in the user guide.

Prior to using the device for the first time, it is recommended that you conduct the following skin test on your forearm:

1) Keep your forearm dry and clean.

2) Turn on the device which will activate level one of the EMS mode, and conduct a three-minute treatment. Gently press device over your skin so it is in full contact. If rash, itch or irritation occurs on the tested skin, please stop the test immediately. If skin feels abnormal, also stop use immediately. 3) If the area is not red and no negative condition or fever appears after 24 hours, then you may keep using the device.

4) It is recommended that you conduct the foregoing skin test prior to using each of three modes prior to usage.

## WARRANTY INFORMATION:

Please use this as your reference for any customer service or warranty needs you may have when the device is purchased within the United States. When purchasing outside the United States, please contact your place of purchase for local warranty information. Our warranty is limited to the terms set out below:

What is Covered in the Warranty:

This lifetime warranty covers only defects in materials and workmanship in your CelluFade BodyBliss when purchased through an authorized partner or directly from LUNA/MOUR. Your warranty is valid for you alone, the original purchaser, and cannot be transferred. To validate your warranty, please keep your original purchase receipt together with these warranty conditions for the duration of the warranty period. If the product was a gift, obtain a copy of the original receipt or have the original purchaser process the warranty claim on your behalf.

What is Not Covered by this Warranty?
This warranty does not cover: wear, tear, or deterioration, damage, failure caused by neglect, abuse or misuse, improper care, impact or insertion of oreign objects, use with a voltage connector, repair or alteration carried out by parties other than LUNA'MOUR or its authorized agents, or loss of device.

What Will LUNA'MOUR do if your device is defective? What Will LUNA'MOUR do if your device is defective?

If you discover a defect in your device covered by the warranty, contact
Customer Care at hello@lunamour.co for assistance. If we confirm a defect
covered by this warranty, LUNA'MOUR will at its sole discretion decide whether
to replace your device with either a new or refurbished device. If we find that
your device is not defective, as determined by LUNA'MOUR in its sole discretion,
we will notify you and return the device to you. All warranty services must be
authorized and performed by LUNA'MOUR or a service provider authorized by
ILINA'MOUR LUNA'MOUR.

WARRANTY LIMITATION & DISCLAIMERS This warranty gives you specific legal rights and is in addition to your statutory rights. Your remedies for the breach of this limited warranty are limited solely to the repair and replacement of your device as set forth in this warranty statement.

To the maximum extent permitted by applicable law, LUNA'MOUR disclaims all implied warranties. In no event shall LUNA'MOUR be liable for any special, incidental, consequential, punitive, or indirect damages based upon breach of warranty, breach of contract, tort (including negligence or strict product liability), or otherwise, even if LUNA'MOUR is advised or aware that such damages may arise. Such damages include but are not limited to: loss of device, cost of any substitute equipment, the claims of third parties, or damage to property.

If the applicable law does not allow limits on warranties, or on the remedies for breach in certain transactions, the limits in this paragraph and the disclaimer or exceptions from warranty may not apply. This limitation does not apply to claims for death or personal injury caused by a defect in the device or any other liability that cannot be excluded or limited by law.

These warranty terms, conditions, and limitations allocate the risks of produc failure between LUNA 'MOUR and you, the purchaser of a LUNA'MOUR branded device. This allocation is recognized by both parties and is reflected in the price of LUNA'MOUR branded device. By using the LUNA'MOUR branded device. you acknowledge that you have read these terms, conditions, and limitations understand them, and are bound by them.

- 9 -

CelluFade BodyBliss

EMS + Vacuum Sculpting+ Red Light Therapy Body & Face Device

HOW TO USE VIDEOS.



Client Care: For information regarding device usage, warranty, and additional information, please contact Client Care at hello@lunamour.co

This product is not intended to diagnose, cure or prevent any disease. LUNA MOUR makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin conditions. A qualified health professional should be consulted with regards to any condition requiring medical attention.

Item Number: II-CFIIII

- 10 -

- 5 -

- 8 -